

## **Before Installing Mini-X Unsigned Drivers in Windows 8**

- Windows 8 by default will not let you install unsigned drivers.
- In order to install an unsigned driver you can follow the steps shown below.
- The steps shown below start Windows 8 in "Disable Driver Signature Enforcement" mode.

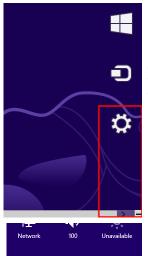

Move the cursor over bottom right corner of the screen. This brings up the **Charm Bar** in Windows 8. Choose the Settings option (icon looks like a Gear).

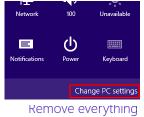

If you want to recycle your PC or s

Advanced startup

Start up from a device or disc (sucl startup settings, or restore Windov

Get started

Restart now

Choose Change PC settings option




Choose **General** on the left hand side and then click **Restart now** button on the right




Click Troubleshoot



**Click Avanced options** 

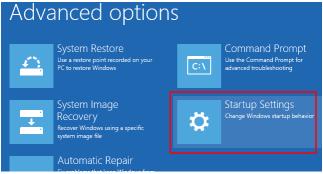

**Click Startup Settings** 

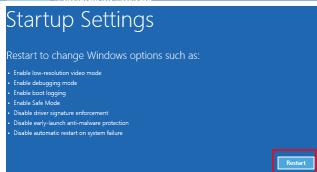

Click the **Restart** button

## Startup Settings

- 1) Enable debugging
- 2) Enable boot logging
- 3) Enable low-resolution video

- 7) Disable driver signature enforcement
  8) Disable early launch anti-malware protection

Press F10 for more options Press Enter to return to your operating system Choose the **Disable Driver Signature Enforcement** and hit the Enter key to start Windows

You are now ready to install the Mini-X drivers. After installing the Mini-X drivers, reboot the computer. Windows 8 will back in "Driver Signature Enforcement" mode.